
GERMAN  
POOL

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# eye 眼部按摩器 massager | me11

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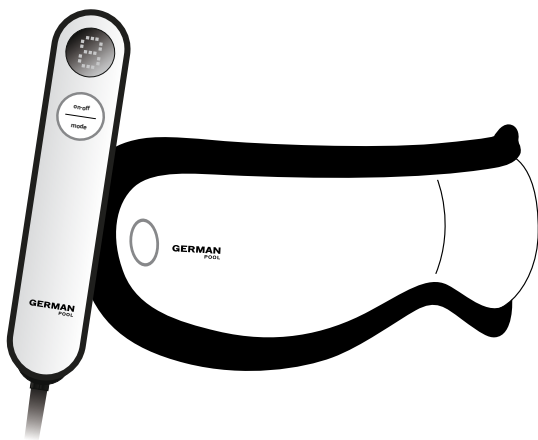---

## IMPORTANT

For safety purpose, you should read the instructions carefully before using the appliance for the first time. Save the instructions for future reference.

The appliance is designed for domestic use and may only be operated in accordance with these instructions for use.

1. Please do not use the eye massager when driving.
2. Please do not use this product in the environment with high temperature or humidity.
3. Make sure contact lenses are removed before using.
4. Wear the device and then switch it on to avoid damaging the air bag due to overpressure.
5. Do not use this product if you suffer from eye diseases, such as glaucoma, cataract and retinal detachment.
6. Do not use this product if you have undergone any eye surgery.
7. Do not use this product if you suffer from strokes, hypertension or hypotension.
8. Sluggish patients, or people with physical or mental or intellectual disability are not allowed to use this massager.
9. Please use the DC 4.5V power adapter provided and pay attention to the safe use of electricity. Battery cannot be put into the battery box if using the external power supply.
10. Please take out the battery if its power is not enough or it will not be used for a long time, so as to avoid damage to the product due to leakage of the battery fluid. Please put waste battery into designated recycling bin to avoid environment pollution. (Alkaline battery is recommended.)
11. This massager is designed for household use, it is not suitable for diagnosis, treatment and other commercial uses.
12. If cleaning is required, wipe it with cotton soaked with alcohol. Do not use water to wash the device.
13. When folding the main unit, never put anything in the center, or touch the center with your finger or any part of the body, to avoid damage or injury.
14. The massager has heating function, users who are insensitive to temperature should pay extra attention.


# INSIDE THE PACKAGE

---

This massager comes with a main unit, a corded remote control, a power adapter and operating instructions

Main Unit

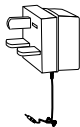Power Adapter

Operating Instructions


Corded Remote Control


# PRODUCT STRUCTURE

---

## Main Unit
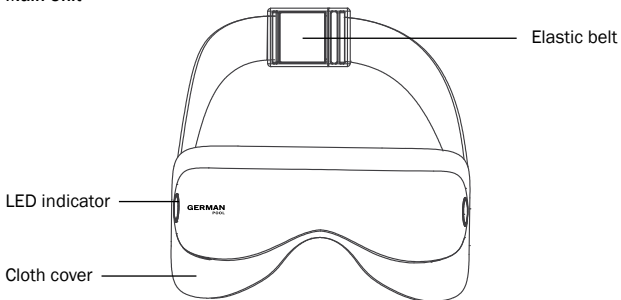
## Corded Remote Control
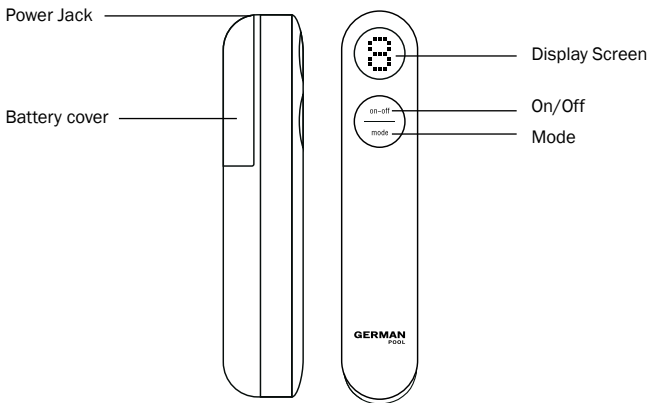
# OPERATING INSTRUCTIONS

---


## Preparation

1. Check if the massager and accessories are intact.
2. Check connection between corded remote control and main unit.
3. Check if there is enough power in the batteries.
4. Check if there is no damage on power adapter.

## Installing batteries

Open the battery cover at the top of the corded remote control, install 3 AA batteries in accordance with the anode and cathode indication. Close the battery cover.

(as in the right figure)
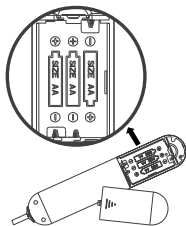
## Adapter

Insert the plug of the power adapter into the power jack of the corded remote control, while the other end is connected to power supply.

## Putting on the massager

Please wear the eye massager with your eyes against the cloth cover of the main unit. Adjust buckles of the elastic belt to wrap the massager around your head comfortably.

## Switching on and off

Press and hold the on-off button for 3 seconds.

## Default

Low Pressure(1) + Low Heat(h) + 15-minute timer

## Manual Mode

Heat Selection -

Press the on-off button to select Low Heat(h), High Heat(H) or Heating off.

Air Pressure Selection -

Press the mode button to select Low Pressure(1), Medium Pressure(2) or High Pressure(3)

# CARE & MAINTENANCE

---

## Power Check

- Check if battery is low or not installed properly. If power adapter is to be used, check the power supply and do not switch on the device if the voltage exceeds the rated range of the device.
- Check the eyes cover and cable regularly. Stop using the device immediately if there is any breaking or rupture. Please send it to German Pool repair center.

## Maintenance

- Take measures to prevent moisture, rats and pollution.
- Switch off the device before cleaning.
- Put the device in safe place after using. Do not put heavy objects on the device.
- Contact qualified technician if there is any fault. Do not disassemble the device yourself.

# TROUBLESHOOTING

---

If there is any problem when you use the device, please follow the instructions below to check. If the problem still exists or the device cannot work normally, please turn off the power supply and contact service center. Do not disassemble or repair the device yourself.

Do not disassemble or repair the device yourself.

Phenomenon	Possible Cause	Solution
Power-on fail.	Whether the air hose is damaged or not connected properly.	Check the condition of the air hose.
	Whether the battery is low.	If yes, change battery. Make sure the battery is installed and fixed in the battery holder.
	Whether the battery is incorrectly installed.	Make sure the +,- end of the battery is installed correctly and fixed in the battery holder.
	Whether the power adapter is connected properly.	Please ensure it is connected properly.
Air Pressure and heat are too low.	Whether the battery is low.	If yes, change battery. Make sure the battery is installed and fixed in the battery holder.

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產品規格及設計如有變更，恕不另行通知。  
Specifications are subject to change without prior notice.

中英文版本如有出入，一概以中文版為準。  
If there is any inconsistency or ambiguity between the English version and the Chinese version, the Chinese version shall prevail.

請填寫以下空格並與發票一併保存，以便維修時核對資料。用戶提供本公司之所有資料必須正確無誤，否則免費保用將會無效。  
Please fill out the form below. This information and the purchase invoice will be required for any repairs. Warranty will be invalid if information provided is found to be inaccurate.

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## Warranty Terms & Conditions

Under normal operations, the product will receive a 1-year full warranty provided by German Pool (effective starting from the date of purchase).

\* This warranty is not valid until Customer's registration information is returned to and received by our Service Centre within 10 days of purchase via one of the following methods:

- 1) Visit our website and register online:  
[www.germanpool.com/warranty](http://www.germanpool.com/warranty).
- 2) Complete all fields on the attached Warranty Card and mail it back to our Customer Service Centre along with a photocopy of the purchase invoice (Warranty Card and invoice copy will not be returned to Customer).

1. Customer who fails to present original purchase invoice will be disqualified for free warranty.
2. Customer should always follow the operating instructions. This warranty does not apply to:
  - labour costs for on-site installation, check-up, repair, replacement of parts and other transportation costs;
  - damages caused by accidents of any kind (including transport and others);
  - operating failures resulted from incorrect voltage, improper operations and unauthorized installations or repairs.
3. This warranty is invalid if:
  - purchase invoice is modified by unauthorized party;
  - product is used for any commercial or industrial applications;
  - product is repaired or modified by unauthorized personnel, or unauthorized parts are installed;
  - serial number is modified, damaged or removed from the product.
4. German Pool will, in its discretion, repair or replace any defective part.
5. This warranty will be void if there is any transfer of ownership from the original purchaser.

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© All rights reserved. Copying, reproducing or using the contents of this manual is not allowed without prior authorization from German Pool, violators will be prosecuted.

[www.germanpool.com](http://www.germanpool.com)